

Product Overview

ProtPure Q Agarose resins are strong anion exchangers based on a 6% highly cross-linked agarose matrix.

They are designed for high-resolution purification of proteins and other biomolecules, offering excellent chemical stability, reproducibility, and scalability for process development and manufacturing.



High Binding Capacity
Consistent & Reproducible



Robust Chemical Stability
Wide pH & CIP Range



Scalable Performance
From Lab to Manufacturing



Agarose-Based
Bioprocess Platform


Comparative Performance Overview

Parameter	Q Agarose Fast Flow	Q Agarose Precise	Q Agarose HR	Typical Market Specifications
BioProcess resin	Yes	Yes	Yes	Yes
Ion Exchanger Type	Strong anion exchanger	Strong anion exchanger	Strong anion exchanger	Strong anion exchanger
Ionic Capacity	0.18–0.25 mmol (Cl ⁻)/mL resin	0.18–0.25 mmol (Cl ⁻)/mL resin	0.18–0.25 mmol (Cl ⁻)/mL resin	0.18–0.25 mmol (Cl ⁻)/mL resin
Ligand	Quaternary Amine	Quaternary Amine	Quaternary Amine	Quaternary Amine
Matrix	6% spherical cross-linked agarose	6% spherical cross-linked agarose	6% spherical cross-linked agarose	6% cross-linked agarose, spherical
Particle Size Range	45–165 µm	45–105 µm	25–55 µm	45–165 µm
Particle Size (d50V)	~90 µm	~65 µm	~35 µm	90 µm
Dynamic Binding Capacity	≈80 mg BSA/mL resin	≈100 mg BSA/mL resin	≈120 mg BSA/mL resin	40–70 mg BSA/mL resin (Fast Flow)
Pressure/Flow Specification	250–450 cm/h, 0.1 MPa, 15 cm bed height	100–300 cm/h, 0.1 MPa, 15 cm bed height	80–150 cm/h, 0.1 MPa, 15 cm bed height	300–600 cm/h, 0.1 MPa, 15 cm bed height
pH Stability (CIP)	2–14	2–14	2–14	2–14
pH Stability (Operational)	2–12	2–12	2–12	2–12
Chemical Stability	Stable in 1.0 M NaOH, 8 M urea, 6 M guanidine hydrochloride, 70% ethanol	Stable in 1.0 M NaOH, 8 M urea, 6 M guanidine hydrochloride, 70% ethanol	Stable in 1.0 M NaOH, 8 M urea, 6 M guanidine hydrochloride, 70% ethanol	Stable in 1.0 M NaOH, 8 M urea, 6 M guanidine hydrochloride, 70% ethanol
Storage Conditions	4–30°C, 20% Ethanol	4–30°C, 20% Ethanol	4–30°C, 20% Ethanol	4–30°C, 20% Ethanol



Applications

- Capture and purification of proteins with basic isoelectric points
- Purification of monoclonal antibodies and antibody fragments
- Viral and plasmid DNA purification
- Separation of host cell proteins and impurities
- Polishing steps in multi-step purification processes
- Process development and scale-up for biopharmaceuticals

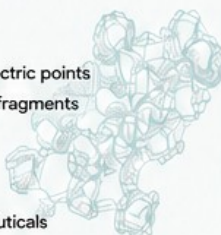

Key Benefits

- ✓ High dynamic binding capacity for efficient capture
- ✓ Excellent resolution and selectivity
- ✓ High flow capacity for faster processing
- ✓ Robust performance across wide pH and cleaning conditions
- ✓ Lot-to-lot consistency for reliable scale-up
- ✓ Designed for process robustness and scalability



Evaluation samples and application support available upon request.
Contact us to discuss your process requirements.

